



INSIGHT
Medical Products Ltd

PRODUCT CATALOGUE

2020

For more information or to order, call

our Sole Distributor in Ireland

Pat Coakley
Medical 

+353 1 2018070
sales@patcoakleymedical.ie



AN INSIGHT INTO INSIGHT

Insight Medical Products Ltd was established by its current directors in May 2000. From our base in the Cotswolds, Insight has become a highly respected supplier of good value, high quality innovative products to the healthcare market.

We are pleased to work with healthcare professionals to design and manufacture innovative products that will help budgets and improve patients well-being.

Insight now has its own in-house manufacturing facility with a fully accredited cleanroom facility for manufacturing and packing. All of our devices comply to the Medical Device Directive and our regulatory systems meet with the requirements of harmonised standards EN ISO 13485 and Directive 93/42/EEC under the supervision of SGS.

We pride ourselves on our reputation for Customer Service and are committed to ensuring a first class service to all of our customers, in the UK and overseas.



CONTENTS

TRACHEOSTOMY CARE

- Tracheostomy Tube Holder
- Trachie-Voice Tracheostomy Speaking Valve
- Tracheostomy Tube Cleaning Swabs
- Advazorb Plus (T) Tracheostomy Dressing
- Advazorb (T) Lite Tracheostomy Dressing
- Advadraw (T) Laminated Tracheostomy Dressing
- Advazorb Fixation Tracheostomy Dressing



TUBE FIXATION

- Endotracheal Tube Holder
- Insil Silicone Tape
- Transducer Holder
- Tube Tidy Clips



LARYNGECTOMY PROTECTORS

- Sofshield Laryngectomy Protectors
- Sofnex Laryngectomy Protectors
- Cascade Shower Protector



LARYNGECTOMY CARE


- Purifoam Laryngectomy Filters
- Actibalm Manuka Honey Lip Balm
- Voice Prosthesis Brush



CONTENTS

AIRWAY MANAGEMENT

- Tracheal Tube Introducers

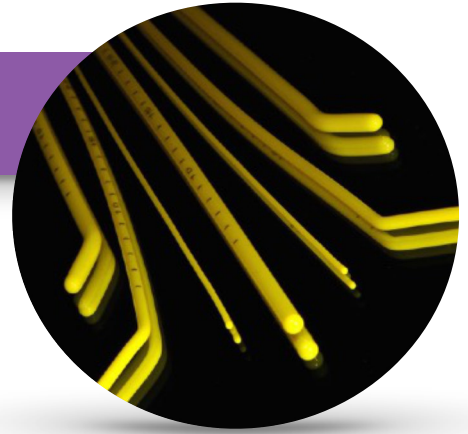

THEATRE CONSUMABLES

- Post Operative Eye Shields
- Insight Bite Guard



MISCELLANEOUS

- Curved Amniotic Hooks
- Female Intermittent Catheters



TRACHEOSTOMY TUBE HOLDER

The Insight Tracheostomy Tube Holder uses the highest quality materials to provide a safe and secure fixation for the tracheostomy tube.

The neck band is manufactured from an ultra soft brushed cotton/foam material that prevents soreness and discomfort to the patient.

A simple to use two-piece design allows quick fitting to the majority of neck sizes, while the integral 'cough strip' enables reflex movement to the patient if required.



- Adult and Paediatric sizes available
- 2-Piece design allows secure fitting for most sizes
- Cough reflex section offers more patient comfort
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
TH/100	Tracheostomy Tube Holder	Adult	289-4814	FTH000	10
TH/150	Tracheostomy Tube Holder	Paediatric	289-4806	FTH092	10

TRACHIE-VOICE TRACHEOSTOMY SPEAKING VALVE

The Trachie-Voice low profile speaking valve fits securely to a 15mm connector tracheostomy tube.


It has a one way mechanism that opens when the patient breathes in, allowing air into the airway via the tracheostomy.

The valve closes on expiration allowing air to follow the normal route of expiration through the mouth and nose, instead of the tracheostomy, permitting speech.

Being 'hands free' allows voicing to take place without the need for finger occlusion, so reducing possible infection and hygiene concerns.



- For non-ventilated patients
- Economical/low cost device
- Lightweight reducing pressure on stoma site
- Reduces hygiene and infection concerns
- Ridged sides enabling quick application and removal
- Supplied complete with storage box
- Easy to clean
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

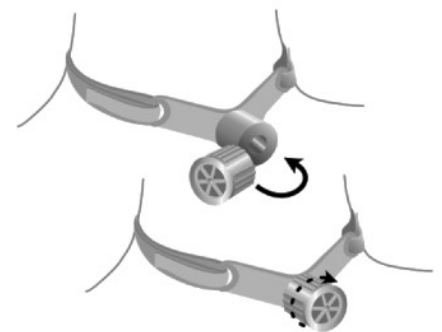

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
TSV/100	Trachie-voice Speaking Valve	One Size	410-7546	FDG1636	12 (6x2)

LATEX
FREE




TRACHEOSTOMY TUBE CLEANING SWAB

The Insight Cleaning Swabs are perfect for cleaning all types of PVC and silicone tracheostomy or laryngectomy tubes.

Available in two head sizes, our Tube Cleaning Swabs are secure and suitable for popular sizes of tubes up to 6mm (small), 9mm (large) and 9mm (Extra Long Handle).

The swabs soft foam heads ensure damaging or scuffing the inside of the tube is avoided, whilst the white plastic handle is flexible for ease of use.

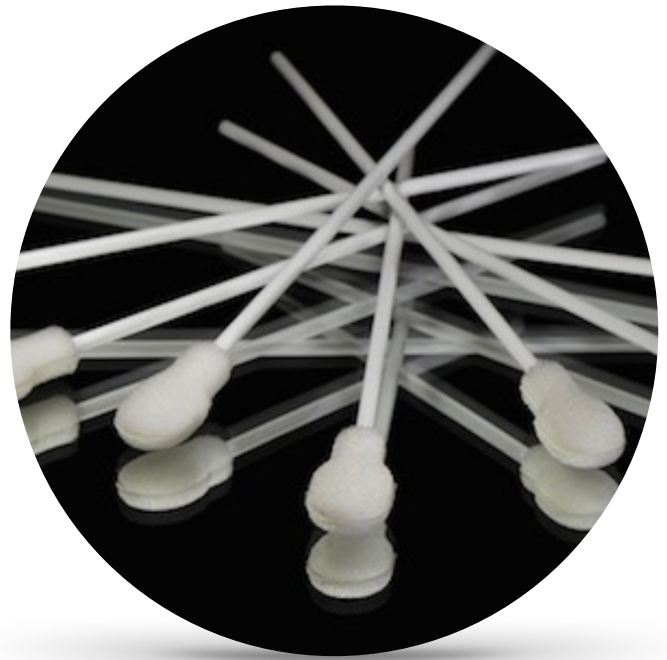

- Use with warm water to remove mucus residue in the tube
- Available in two sizes sizes suitable for popular sizes of tube
- Available on FP10 Prescription
- Single patient use only
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
TCS/900	Tube Cleaning Swabs	Large	324-8689	FDH420	100 (10x10)
TCS/950	Tube Cleaning Swabs	Small	324-8671	FDH421	100 (10x10)
TCS/1000	Tube Cleaning Swabs	Extra Long	N/A	FDH611	100 (10x10)

IMP CODE	HEAD WIDTH	HEAD LENGTH	HANDLE LENGTH	MAXIMUM TUBE SIZE
TCS/900	13mm	26mm	155mm	9mm
TCS/950	8mm	22mm	155mm	6mm
TCS/1000	9mm	25mm	245mm	9mm

ADVAZORB PLUS (T) TRACHEOSTOMY DRESSING



Advazorb Plus (T) offers a good standard in hydrophilic foam tracheostomy dressings with high absorbency ability.

Manufactured in lightweight, soft and conforming low-adherent foam with high fluid absorbance and retention characteristics.

Manufactured with a film backing to ensure improved strength and absorbance for medium/high exuding wounds.

Supplied sterile and individually wrapped, ready for use.

Available in one convenient size to suit most intubation procedures, and to offer comfort to the patient.

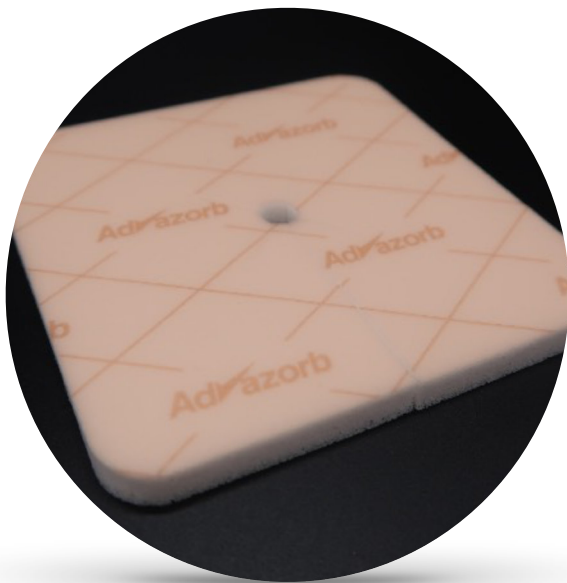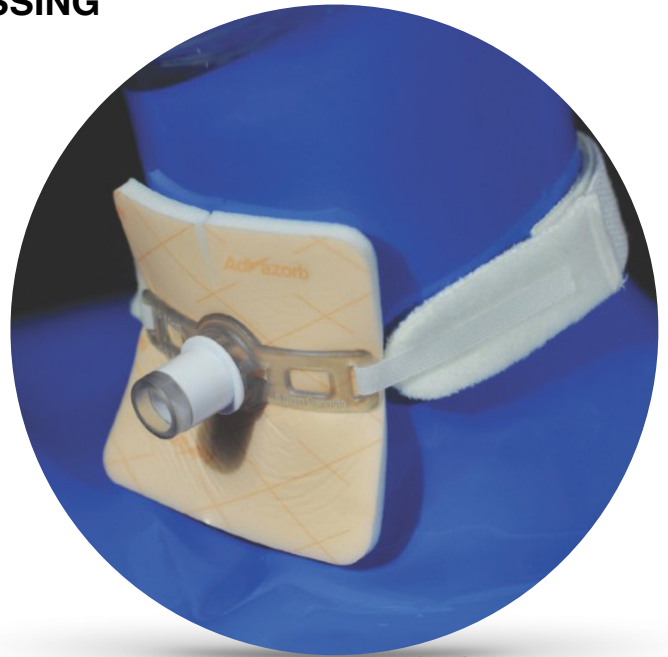

- Cushioning effect on wound
- Low adherence -Reduced trauma on removal
- Highly absorbent 5mm foam
- Excellent patient comfort
- Rapid uptake of exudates
- Can be cut easily to size
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/3900	Advazorb Plus Dressing	80 x 80mm	395-1126	FDH419	20



ADVAZORB (T) LITE TRACHEOSTOMY DRESSING

Advazorb (T) Lite is the definitive lightweight hydrophilic slim foam tracheostomy dressing.

Advazorb (T) offers remarkable standards of conformability, softness, lightness compared to gauze dressings.

Perfect for the treatment of light/medium exuding wounds as moisture is absorbed and retained in the hydrophilic cellular structure.

Supplied sterile and individually wrapped, ready for use. Available in one convenient size to suit most intubation procedures, and to offer comfort to the patient.

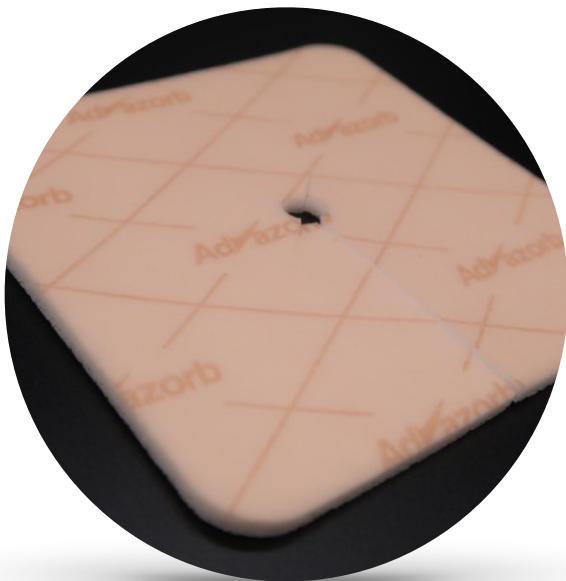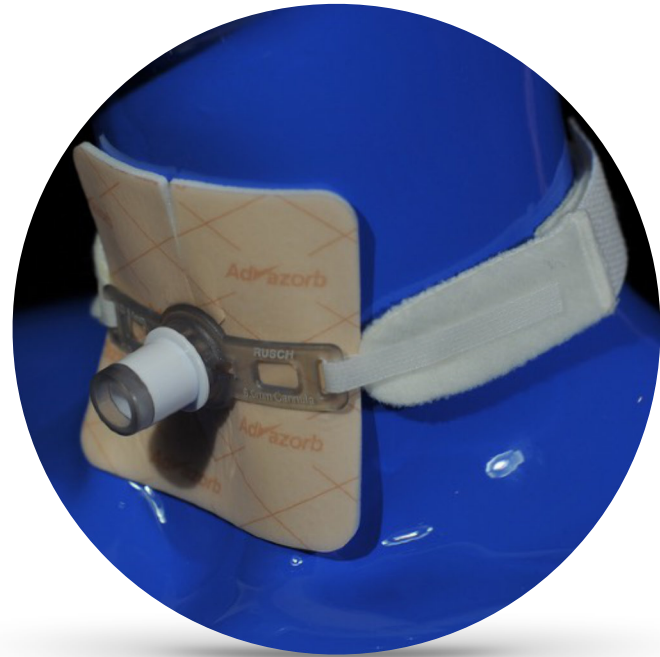

- Cushioning effect on wound
- Low adherence -Reduced trauma on removal
- Highly absorbent 2.5mm foam
- Excellent patient comfort
- Rapid uptake of exudates
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/4212	Advazorb Lite Dressing	80 x 80mm	374-0230	ELY392	20

ADVADRAW (T) LAMINATED TRACHEOSTOMY DRESSING

Advadraw 'T' is a tri-laminate tracheostomy dressing which has a low adherent perforated wound film contact layer, an absorbent centre action and a strike-through resistant layer to seal the dressing.

Advadraw 'T' is comprised of an enhanced lighter fibrous material, with tests showing improved absorbency rates when compared to other manufacturers.

The slim contours of Advadraw 'T' avoids elevating the tracheostomy tube away from the neck, while the curved edges ensure maximum comfort to the patient.

Available in one convenient size to suit most intubation procedures, can be easily trimmed if needed.



- Cushioning effect on wound
- Low adherence -Reduced trauma on removal
- Conformable patient protection
- Suitable for low/medium levels of exudate
- Rapid uptake of exudates
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/4416	Advadraw (T) Dressing	80 x 80mm	409-4970	FDH615	20


ADVAZORB FIXATION TRACHEOSTOMY DRESSING

Advazorb Fixation is a post tracheostomy decannulation dressing fitted with a central location button that prevents unwanted escape of air through the wound site with gentle pressure applied during occasions such as talking, coughing and sneezing.

Advazorb Fixation comprises of a hydrophilic foam dressing with a unique perforated silicone adhesive wound contact layer. It is an atraumatic foam dressing which does not adhere to moist wounds.

The silicone helps maintain a moist wound healing environment and makes it ideal as an adherent wound contact dressing. Soft silicone dressings do not damage newly formed granulating tissue and do not lose their adhesion after initial application.

The silicone layer prevents the need for additional dressing or fixation. Silicone adhesive combines security with pain-free dressing changes. Waterproof backing enables the patient to shower without the need to replace or cover the dressing.

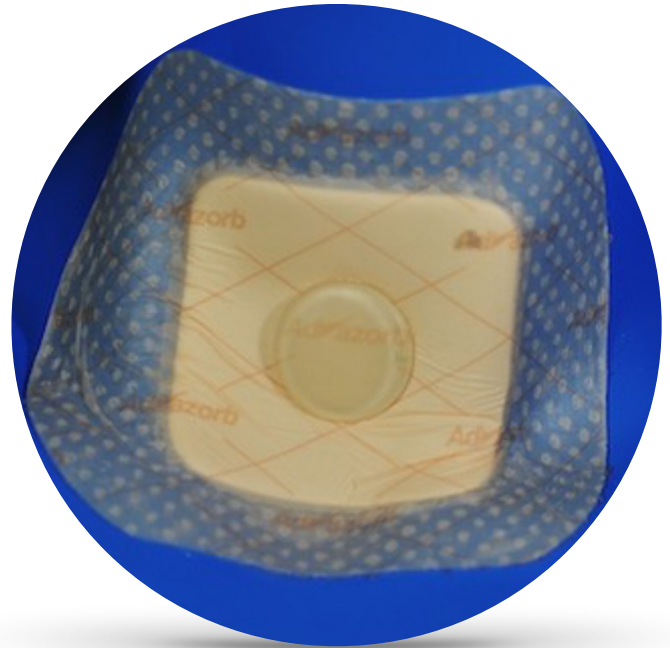

- Supplied sterile and individually wrapped, ready for use
- Location button enables accurate identification of tracheostomy site
- Maintains a moist environment and manages exudates
- Allows atraumatic dressing change
- Provides a barrier for bacteria
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/4341	Advazorb Fixation	12cm x 10cm	N/A	FDH605	10

ENDOTRACHEAL TUBE HOLDER

The Endotracheal Tube Holder Set has been commended by professionals on the ease of use, security and patient comfort.

Because of the soft materials, the product has been proven to greatly reduce the incidents of sores to the patients mouth, while maintaining a secure fitting to the ultra soft neckpiece.

The neckpiece is fitted with an adjustable laminated velcro tab which can be repositioned to size.

Each set is supplied with 3 holders to enable the user to replace soiled items, therefore improving value over alternative methods.



DETAILED INSTRUCTION SHEET AND PRODUCT TRAINING AVAILABLE ON REQUEST



- Fully adjustable to fit most neck sizes
- Two sizes available
- Spare white neckbands available to order
- Laminated hook/loop tabs offer total security of fitting
- Ultra soft materials prevent abrasions to the patient
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	NHS CODE	PACK QTY
TFP/200	Endotracheal Tube Holder Set	Up to 20" neck	FTH005	10
TFP/201	Endotracheal Tube Holder Set Extra Long	Up to 24" neck	FTH025	10
TFP/222	Endotracheal Tube Holder Neck Band Only	Up to 20" neck	FTH130	10
TFP/225	Endotracheal Tube Holder Neck Band Only Extra Long	Up to 24" neck	FTH131	10



LATEX
FREE



INSIL SILICONE TAPE

This versatile, breathable tape is backed with a silicone adhesive. When tape is removed it is extremely gentle on the skin, however it can be re-applied without losing adhesive qualities.

Insil multi-use tape can be used to fix most medical devices such as tubes, IV cannulae, eye shields, dressings and drains.

Can be used in theatres for securing closure of eyelids during operations, and allowing painless removal.

Insil is available in 2cm and 4cm widths with fine perforations every 5cm to allow separation without the need for scissors. Supplied in pocket-sized dispenser boxes for ease of use.



- Available on FP10 Prescription
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
INSIL/200	Insil Silicone Tape	20mm x 3m	342-4132	IMP Direct	12 (12x1)
INSIL/400	Insil Silicone Tape	40mm x 1.5m	342-4140	IMP Direct	12 (12x1)

LATEX
FREE





TRANSDUCER HOLDER

A simple and effective method for ensuring that the transducer remains in the correct position either within the department or while the patient is being moved and replaces the need for adhesive tape application.

The holder ensures stable positioning of the transducer for extended periods of time whilst providing patient comfort.

Made from soft material and secured to the arm with fully adjustable hook and loop fastener to fit most patients, eliminating the difficulties in taping the transducer to the patient or repositioning the transfusion stand due to patient movement, also saving nursing time.

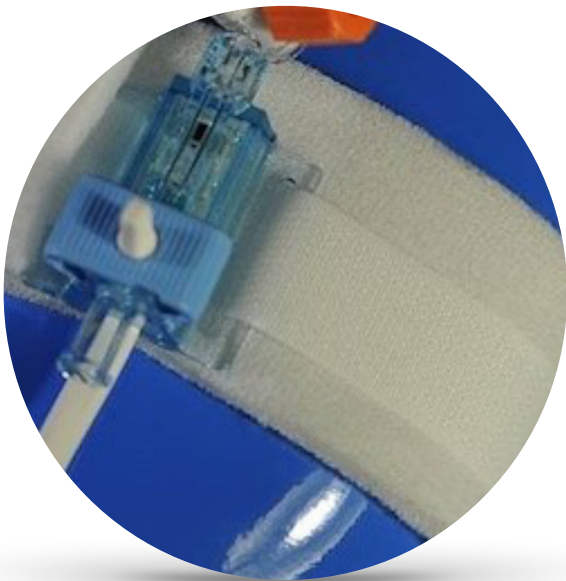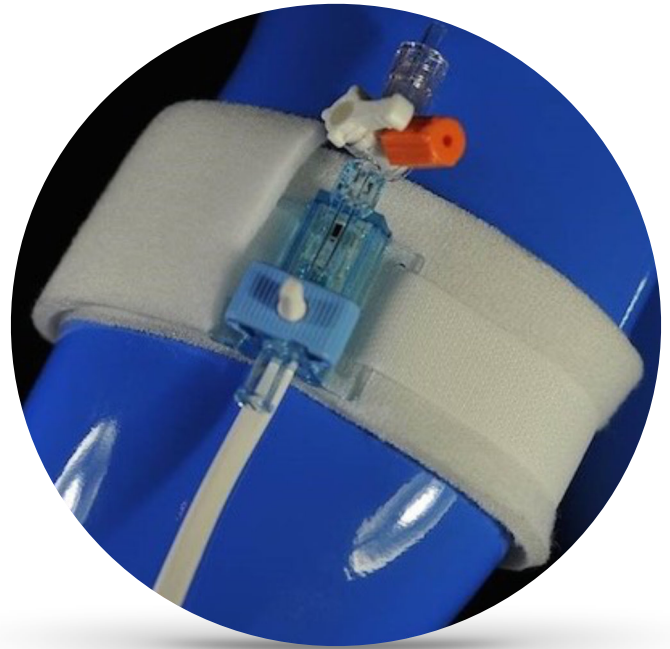

- Designed in partnership with the ICU at the Royal Hallamshire Hospital
- For use in the department or while patient is being transferred
- No need for any use of adhesive tape on the patient
- Quality soft foam/cotton material for patient comfort
- Easy to apply
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
TFP/400	Transducer Holder	Standard	N/A	IMP Direct	10

LATEX
FREE



TUBE TIDY CLIPS

The simplest product on the market for keeping tubes and lines arranged and sorted during patient monitoring.

The product can clip to blankets, sheets, pillows and opens to allow tubing to run through.

Simply lay tubes required over hook/loop and seal to hold tubes in place.

Organises cords and drapes, useful in all specialist units and wards.

Can retain drapes and linens without the need of pins or fixers, thus eliminating the chance of personal injury



- Available in white, blue, red and yellow for tube identification
- User is able to write on the product for identity markings
- Simply open the hook/loop material, insert tubing and close
- Quality soft foam/cotton material for patient comfort
- Easy to apply
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	COLOUR	DRUG TARIFF	NHS CODE	PACK QTY
TFP/800	Tube Tidy Clips	White	N/A	IMP Direct	10
TFP/850.B	Tube Tidy Clips	Blue	N/A	IMP Direct	20
TFP/850.R	Tube Tidy Clips	Red	N/A	IMP Direct	20
TFP/850.Y	Tube Tidy Clips	Yellow	N/A	IMP Direct	20

LATEX
FREE



SOFSHIELD LARYNGECTOMY & TRACHEOSTOMY PROTECTORS



Sofshield Laryngectomy & Tracheostomy Protectors can be washed a minimum of 10 times, while still retaining filtration efficiency.

This offers a significant cost saving and longer product use over alternative protectors that use foam inserts, while also eliminating the need for costly additional cosmetic covers.

All Sofshield Protectors have an antibacterial coating to help prevent infections.

Sofshield is available in five colours and three sizes to suit all Laryngectomy and Tracheostomy patients.

Sofshield uses a silicone toggle system that eliminates the need to tie cords, therefore allowing patients with restrictive movement an easier fitting process.

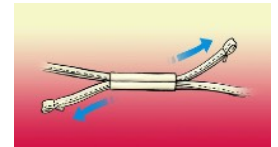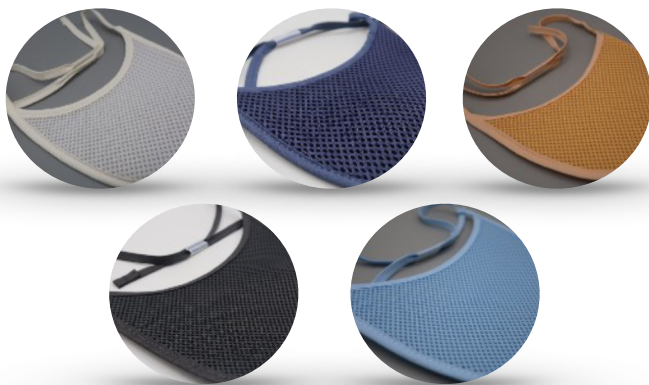

- Antibacterial coating to prevent infection
- Consistent filtration efficiency after being washed
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
AS/3669	Sofshield -White	Large	291-9231	FTH006	10
AS/3764	Sofshield -Blue	Large	305-5373	FTH008	10
AS/3765	Sofshield -Beige	Large	305-5381	FTH007	10
AS/3698	Sofshield -White	Small	291-9249	FTH009	10
AS/3777	Sofshield -Blue	Small	305-5399	FTH100	10
AS/3778	Sofshield -Beige	Small	305-5407	FTH010	10
AS/4051	Sofshield -Black	Medium	377-1722	FSM2621	10
AS/4056	Sofshield -Light Blue	Medium	409-1526	FTH129	10


SOFNEX LARYNGECTOMY PROTECTOR

Sofnex Protectors can be washed a minimum of 10 times, offering a significant cost saving and longer product use over alternative protectors that still use foam inserts whilst still retaining filtration efficiency.

Sofnex is supplied with a velcro fastening feature that allows fixation to the enclosed soft neck straps.

All Sofnex Protectors have an antibacterial coating to help prevent infections and available in three colours with 2 straps per box.

The appearance of Sofnex eliminates the need for costly additional cosmetic covers.

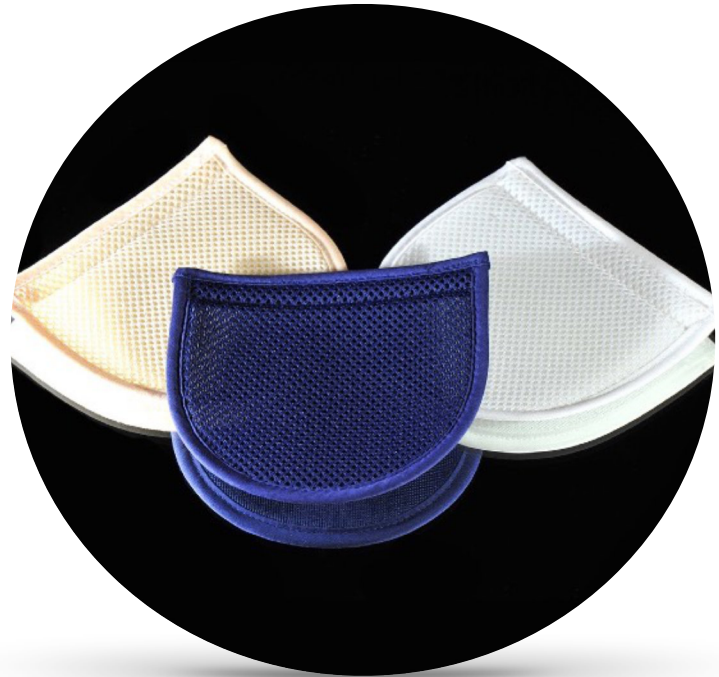

- Antibacterial coating to prevent infection
- Consistent filtration efficiency after being washed
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	COLOUR	DRUG TARIFF	NHS CODE	PACK QTY
CR/3856	Sofnex Protector	White	319-0402	IMP Direct	10
CR/3859	Sofnex Protector	Blue	319-0410	IMP Direct	10
CR/3860	Sofnex Protector	Beige	319-0428	IMP Direct	10

LATEX
FREE



CASCADE SHOWER PROTECTOR

Cascade is a soft waterproof PVC protector that provides a barrier to prevent shower water from entering the stoma thereby improving the patient's confidence to shower routinely.

The neckband is soft and has elastic to form a comfortable yet firm seal around the neck to prevent water seeping below the protector. The backing is a knitted material that will absorb moisture from surrounding wet skin and provide a comfortable warm feel, unlike other products that may initially feel cold to the skin.

From a safety perspective, the material backing also prevents adhesion to the skin or stoma button when bib is wet, as adhesion from unlined protectors can reduce the ability for a patient to breath properly.

For added comfort and breathability the Cascade Protector can be worn over a Sofshield Protector or alternative filter, however, this is not essential due to the backing material. The protector is an aesthetically pleasing white in colour and supplied in pairs to allow



- Shower proof
- Lay flat contour
- Washable
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked


IMP CODE	DESCRIPTION	COLOUR	DRUG TARIFF	NHS CODE	PACK QTY
AS/3835	Cascade Shower Protector	White	314-2601	FTH028	10 (5x2)

PURIFOAM LARYNGECTOMY FILTERS

Purifoam Filters use a high quality medical grade foam that protects the stoma from dirt and dust entering the airway, while warming and humidifying the incoming air.

Purifoam is adhered to the skin using a skin friendly adhesive strip that has rounded edges for added comfort.

Manufactured from beige coloured foam to give a natural aesthetic look. Purifoam has excellent moisture retaining qualities, ensuring good levels of absorbency to give exceptional humidification to the airway for breathing reassurance.

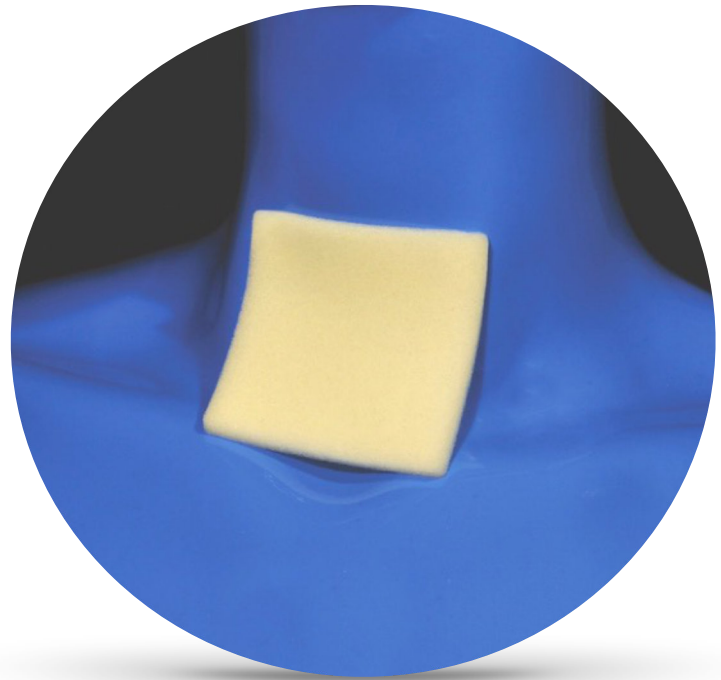

- Change daily or more frequently if required
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/3853	Purifoam Filters	6.5cm x 6.5cm	319-0386	FTH029	30

ACTIBALM MANUKA HONEY LIP BALM

Actibalm provides soothing relief to dry, cracked lips.

Actibalm is a topical ointment made up from only two ingredients; white pharmaceutical grade petroleum jelly and pure Activon Manuka Honey.

Actibalm creates a protective barrier, shielding the skin from exposure to external agents that may cause drying and cracking. This barrier locks in moisture, aiding in the restoration of dry skin.

The product is not for use on broken skin, mucosal surfaces or to be invasive. Do not use the product if allergic to bee venom, bee products or petrolatum.



- Soothing relief
- Moisture protective barrier
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
CR/3955	Actibalm Manuka Honey Lip Balm	10g	347-1224	IMP Direct	20 (20x1)

LATEX
FREE



VOICE PROSTHESIS BRUSHES

Voice Prosthesis Brushes are manufactured with soft nylon bristles that prevent damaging the silicone inner surface of the prosthesis.

Brushes are tapered to provide effective cleaning of mucus and secretion build-up in the inner lumen of the voice prosthesis.

The brush is fitted with a silicone safety sleeve to prevent over-insertion whilst in use.



- Silicone safety sleeve
- Polycarbonate handle for ease of use
- Available on FP10 Prescription
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	DRUG TARIFF	NHS CODE	PACK QTY
PB/50	Voice Prosthesis Brush	Standard	229-4395	IMP Direct	10


LATEX
FREE



TRACHEAL TUBE INTRODUCERS

Tracheal Tube Introducers (Bougies) are designed to help insert tracheal tubes during difficult intubation.

All Introducers are single use to reduce chances of cross-infection and are produced in straight or angled format to suit all procedures.

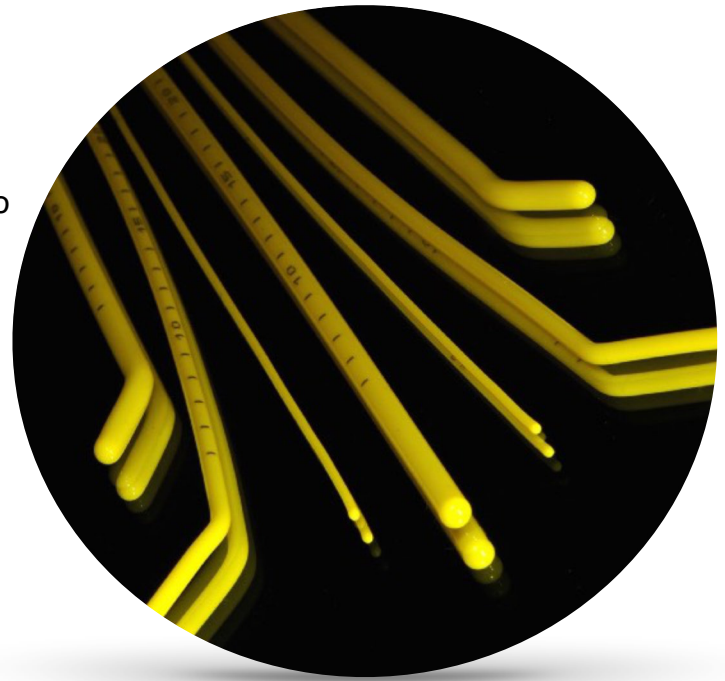

- Sheen surface gives less tracheal tube resistance
- Graduated markings to enhance entry and exit management
- Barium tip to assist X-Ray detection
- Available with straight or angled tip
- Saves costs in cleaning and product deterioration
- Single patient use only
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	LENGTH	NHS CODE	PACK QTY
SB/05/ST	Sterile Bougie -Straight	5ch	470mm	FTH031	10
SB/10/600/S	Sterile Bougie -Straight	10ch	600mm	FTH032	10
SB/10/800/S	Sterile Bougie -Straight	10ch	800mm	FTH033	10
SB/15/600/S	Sterile Bougie -Straight	15ch	600mm	FTH034	10
SB/15/800/S	Sterile Bougie -Straight	15ch	800mm	FTH035	10
SB/10/600/A	Sterile Bougie -Angled	10ch	600mm	FTH036	10
SB/10/800/A	Sterile Bougie -Angled	10ch	800mm	FTH037	10
SB15/600/A	Sterile Bougie -Angled	15ch	600mm	FTH040	10
SB/15/800/A	Sterile Bougie -Angled	15ch	800mm	FTH041	10
SBV/14/750/OX	Angled, Vented with Oxy Con	14ch	750mm	FTH042	10

INSIGHT BITE/TEETH GUARD

Insight Bite Guards are a protective shield used to protect teeth against accidental damage from contact with rigid instruments during surgery.

Enhanced design to provide safer retention when fitted. Manufactured from medical grade TPE for pliability and patient comfort.

Coloured bright yellow to allow operator visibility.




- Shaped for optimum fit to patient
- Radio-opaque, autoclavable on the 134 degree cycle up to 3.5 minutes
- Can be supplied sterile or unsterile in easy peel blister pack
- Single patient use
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	SIZE	COLOUR	NHS CODE	PACK QTY
BG/100	Unsterile Bite Guard	Adult	Yellow	FTH124	20
BG/200	Sterile Bite Guard	Adult	Yellow	FTH125	20

POST OPERATIVE EYE SHIELDS

Eye Shields are used after surgery to protect the eye and to hold dressings in place while holes in the shield provide ventilation.

Supplied bulk for trade suppliers of Theatre Procedure Packs. We also supply individually sterile blister packed for hospitals and can use customers own artwork if required.

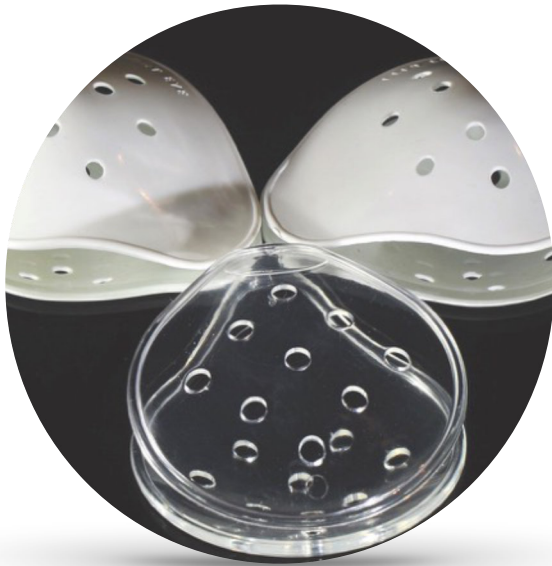

- Universal fitting available for dual fitting
- Available in Clear, Blue Tint or Opaque White
- Left and Right fitting only available in Opaque White
- Single patient use
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE	DESCRIPTION	PACKAGING	COLOUR	PACK QTY
ES/500/C	Universal Eye Shield	Bulk	Clear	2000
ES/500/W	Universal Eye Shield	Bulk	White	2000
ES/500/B	Universal Eye Shield	Bulk	Blue Tint	2000
ES/100/C	Universal Eye Shield	Sterile	Clear	100
ES/100/W	Universal Eye Shield	Sterile	White	100
ES/100/B	Universal Eye Shield	Sterile	Blue Tint	100
ES/500/WR	Right Eye Shield	Bulk	White	2000
ES/500/WL	Left Eye Shield	Bulk	White	2000
ES/100/WR	Right Eye Shield	Sterile	White	100
ES/100/WL	Left Eye Shield	Sterile	White	100

LATEX
FREE




CURVED AMNIOTIC HOOKS

Curve-Hook Anatomically Shaped Amniotic Membrane Perforator.

The Curve-Hook is anatomically shaped to make the hook more versatile and easier to use. The Curve-Hook has been designed to make the membrane easily accessible.

The design ensures optimum efficiency and the prevention of injury to the patient when used.

The shaped rounded tip and protected hook simply engages and ruptures the membrane, saving the need to prod or puncture, therefore reducing any chance of injury to the foetus.

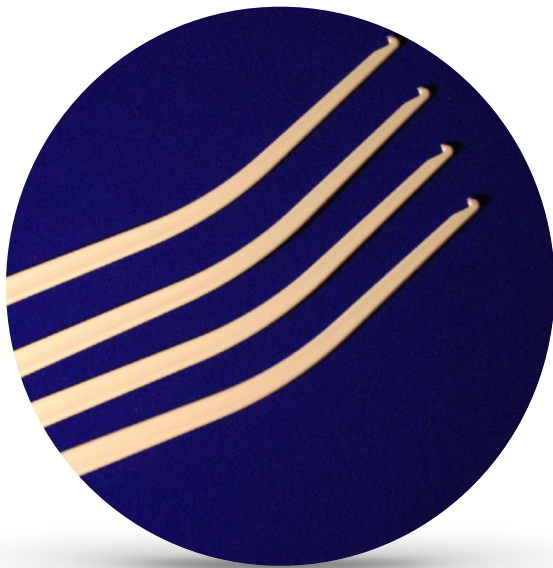

- Individually supplied sterile in peel pouch
- Single patient use
- UK Manufactured
- Latex-Free
- CE Marked

IMP CODE

DESCRIPTION

SIZE

NHS CODE

PACK QTY

KH/100

Curve-Hook

Standard

IMP Direct

100

LATEX
FREE



HANDHELD INTERMITTENT CATHETERS

The Handheld Intermittent Catheter is fully adjustable for comfort and ease of use for self catheterisation.

Movable handle to ascertain the ideal fitting position.

The catheter has a nelson tip and made from PVC to minimise allergic reactions.



- Individually supplied sterile in peel pouches
- Single patient use
- Latex-Free
- CE Marked

IMP CODE	HANDLE SIZE	CATHETER SIZE	NHS CODE	PACK QTY
INS/100/10	115mm Length	180mm x 10ch	IMP Direct	100
INS/100/12	115mm Length	180mm x 12ch	IMP Direct	100
INS/100/14	115mm Length	180mm x 14ch	IMP Direct	100

LATEX
FREE





Pat Coakley 
Medical

For more information or to order, call

our Sole Distributor in Ireland

Pat Coakley Medical Ltd

Units E3 E4
Network Enterprise Park
Kilcoole A63 H244
Co Wicklow
IRELAND

+353 1 2018070
sales@patcoakleymedical.ie

